

Ordering Information

Titanium Ligating Clips

Product No.	Clip Size	Clips/Cartridge	Cartridges/Box
20130	Small, Narrow	6	30
20425	Small, Narrow	24	25
30130	Small, Wide	6	30
30425	Small, Wide	24	25
40130	Medium	6	30
40425	Medium	24	25
50120	Medium Large	6	20
50415	Medium Large	24	15
60120	Large	6	20
60415	Large	24	15

Open Surgical Ligating Appliers

Product No.	Clip Size	Length	Jaw Type
70262	Small, Narrow	6" 15 cm	Curved
70282	Small, Narrow	8" 20 cm	Curved
70285	Small, Narrow	8" 20 cm	55°
70212	Small, Narrow	11" 28 cm	Curved
70362	Small, Wide	6" 15 cm	Curved
70382	Small, Wide	8" 20 cm	Curved
70385	Small, Wide	8" 20 cm	55°
70312	Small, Wide	11" 28 cm	Curved
70462	Medium	6" 15 cm	Curved
70482	Medium	8" 20 cm	Curved
70485	Medium	8" 20 cm	55°
70412	Medium	11" 28 cm	Curved
70415	Medium	11" 28 cm	55°
70417	Medium	11" 28 cm	70°
70582	Medium Large	8" 20 cm	Curved
70585	Medium Large	8" 20 cm	55°
70512	Medium Large	11" 28 cm	Curved
70515	Medium Large	11" 28 cm	55°
70682	Large	8" 20 cm	Curved
70685	Large	8" 20 cm	55°
70612	Large	11" 28 cm	Curved
70615	Large	11" 28 cm	55°
70617	Large	11" 28 cm	70°
70060	Remover-Universal	6" 15 cm	Straight
70010	Remover-Universal	11" 28 cm	Curved


Endoscopic Appliers

Product No.	Clip Size	Length	Diameter	Jaw Type
80502	Medium Large	12" 33 cm	10/11 mm	Curved
90502	Medium Large	17" 45 cm	10/11 mm	Curved
80602	Large	12" 33 cm	12 mm	Curved
90602	Large	17" 45 cm	12 mm	Curved



MEDCHEM
MED-CHEM SURGICAL LTD
 Unit 4 / 118 Asquith Ave
 St Lukes, Auckland 1025
 PO Box 37400
 Parnell, Auckland 1151
 P: 0800MEDCHEM (0800 633 2436)
 (09) 815 2975
 F: (09) 815 2974
www.medchemsurgical.co.nz

Vesocclude Ligation System

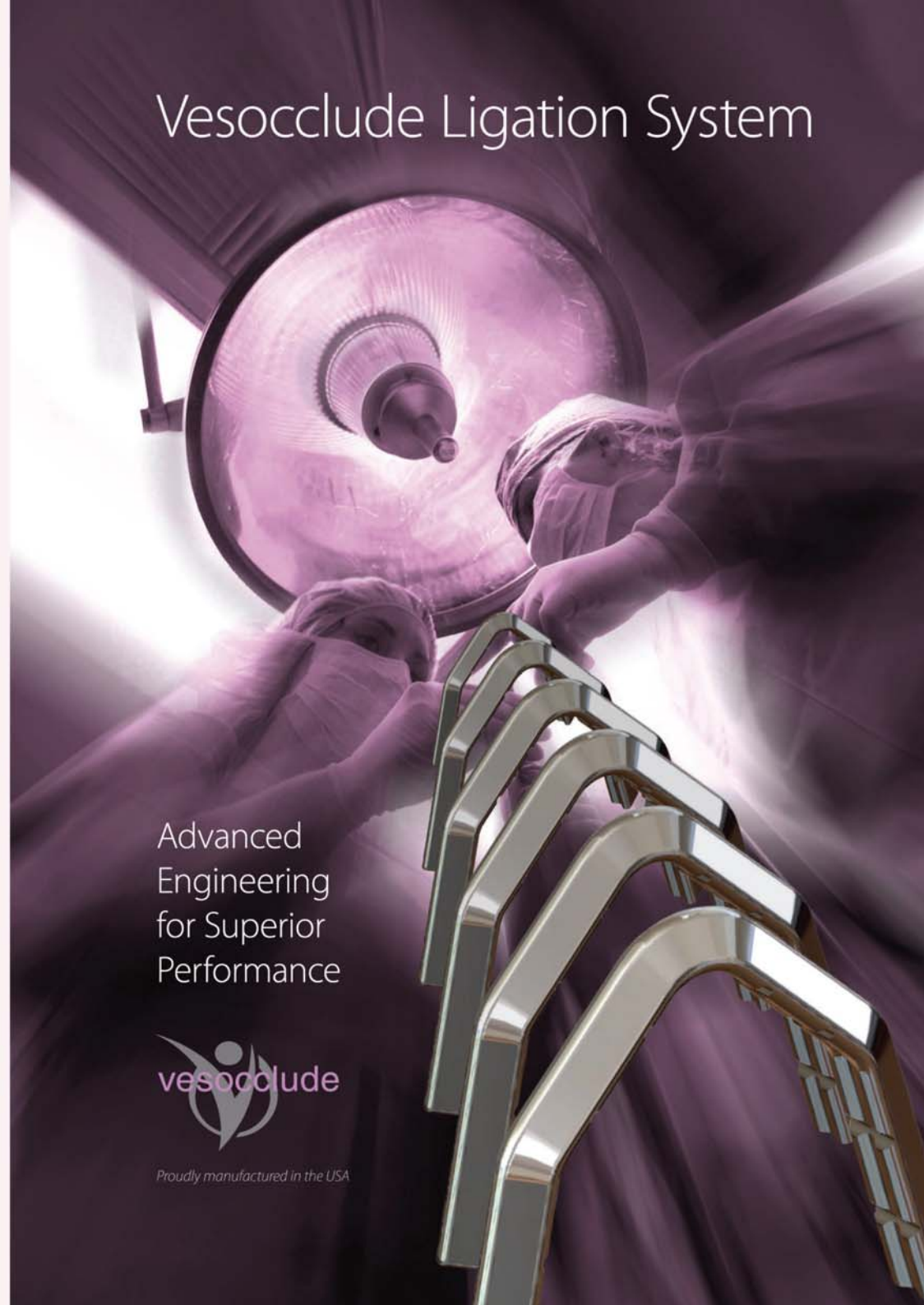

Advanced
 Engineering
 for Superior
 Performance



Proudly manufactured in the USA

Advanced Cartridge System

Loading is easy with our push and click system. You get feedback when a clip is loaded.

Fingers provide a guide for ease of loading clip into applicator. After loading, fingers remain in down position to signify clip has been loaded.*


Adhesive tape on base allows adherence to any sterile surface

Tab lock outs and large base with deep side grooves for ease of handling

Cartridge system contains barium to improve detection

Cartridges are color-coded to match applicators, ensuring proper loading and clip placement

Available in 6 and 24 packages for economical usage

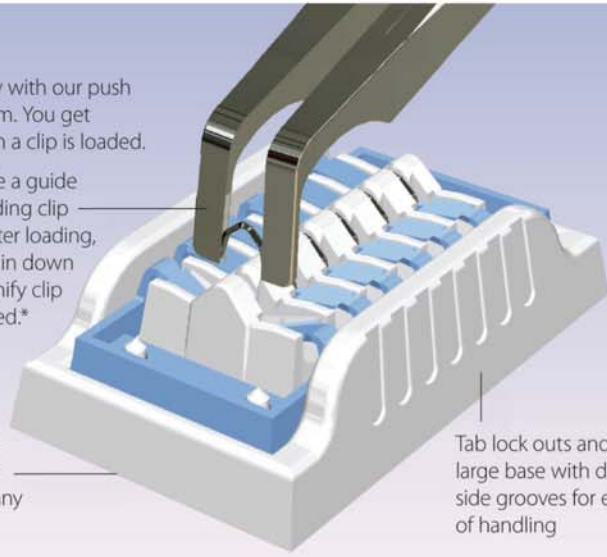

Superior Clip Design

Chevron shape engulfs tissue with precise tip-to-tip closure

Heart-shaped wire ensures firm grip on vessels



Heart shaped cross section maximizes surface-to-surface contact with applicator jaw, practically eliminating clip fallout



Curved




55°



70°

Versatile Applicators for Open and Endoscopic Procedures

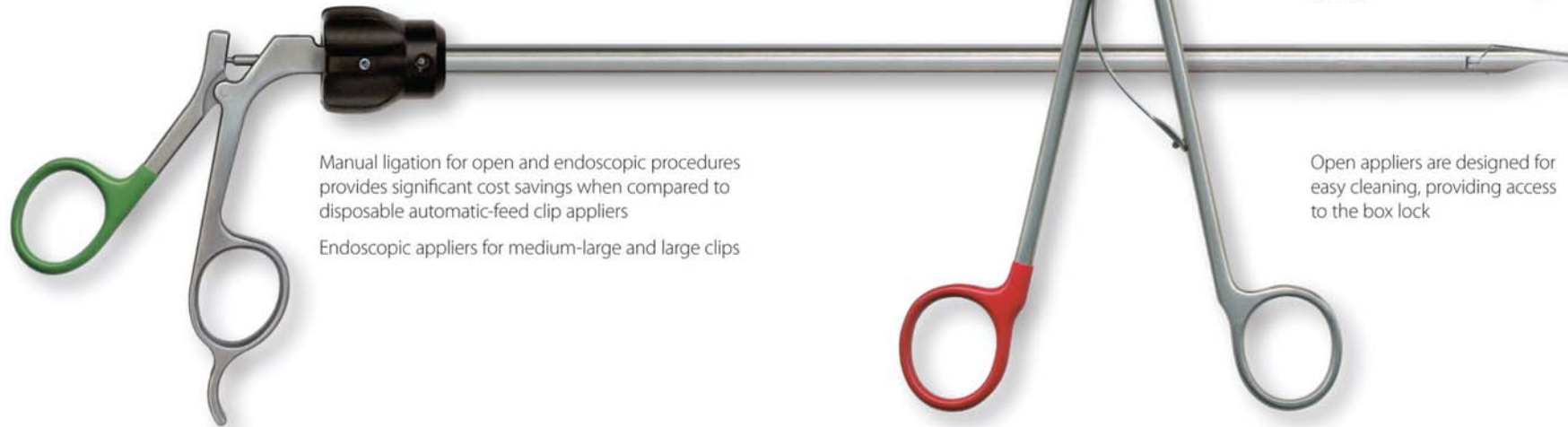

Manual ligation for open and endoscopic procedures provides significant cost savings when compared to disposable automatic-feed clip applicators

Endoscopic applicators for medium-large and large clips

Open applicators are designed for easy cleaning, providing access to the box lock

